

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
7 February 2002 (07.02.2002)

PCT

(10) International Publication Number
WO 02/10771 A1

(51) International Patent Classification⁷: **G01N 33/86,**
33/543

(21) International Application Number: PCT/US01/24132

(22) International Filing Date: 1 August 2001 (01.08.2001)

(25) Filing Language: English

(26) Publication Language: English

(30) Priority Data:
60/222,204 1 August 2000 (01.08.2000) US
60/299,129 18 June 2001 (18.06.2001) US

(71) Applicant (for all designated States except US): **CARE-
POINT DIAGNOSTICS, INC.** [US/US]; 10180 Viking
Drive, Eden Park, MN 55344 (US).

(72) Inventors; and

(75) Inventors/Applicants (for US only): **HENRY, Michael,**
R. [US/US]; 535 Mitchell Drive, Grayslake, IL 60030
(US). **LEAHY, David, C.** [US/US]; 820 Cambridge Drive,
Grayslake, IL 60030 (US). **SHEBUSKI, Ronald, J.**
[US/US]; 918 East Shore Road, Bergland, MI 49910 (US).

(74) Agent: **CHAE, Mark, L.**; McDonnell Boehnen Hulbert
& Berghoff, Suite 3200, 300 South Wacker Drive, Chicago,
IL 60606 (US).

(81) Designated States (*national*): AE, AG, AL, AM, AT, AU,
AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU,
CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM,
HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK,
LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX,
MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL,
TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, ZA, ZW.

(84) Designated States (*regional*): ARIPO patent (GH, GM,
KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian
patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European
patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE,
IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF,
CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD,
TG).

Published:

- with international search report
- before the expiration of the time limit for amending the
claims and to be republished in the event of receipt of
amendments

For two-letter codes and other abbreviations, refer to the "Guid-
ance Notes on Codes and Abbreviations" appearing at the begin-
ning of each regular issue of the PCT Gazette.

(54) Title: ANALYSIS OF BIOLOGICAL SAMPLES FOR PLATELET ACTIVATION OR COAGULATION ACTIVATION
MARKERS USING MICROPARTICULES

(57) Abstract: The invention disclosed herein comprises methods for analyzing a biological sample, such as undiluted or diluted whole blood, as well as fractions thereof, for the presence or absence and/or the concentration of disease-specific and/or other medical condition-specific markers. Such markers may include platelet activation and coagulation activation markers. The methods may comprise combining the biological sample with a coated solid phase and analyzing for the presence or absence and/or the concentration of the markers. The analysis may be performed either before or after separation of the solid phase from the biological sample. The analysis may be performed on the combined components or on any of the separation components. A preferred solid phase may be paramagnetic microparticles coated with antibodies or proteins specific for platelet activation and/or coagulation activation markers.

BEST AVAILABLE COPY

ANALYSIS OF BIOLOGICAL SAMPLES FOR PLATELET ACTIVATION OR COAGULATION ACTIVATION
MARKERS USING MICROPARTICULES**CROSS REFERENCE TO RELATED APPLICATIONS**

5

This application claims priority from U.S. Provisional Patent Application No. 60/222,204, filed 01 August 2000 (Atty. Docket No. MBHB00-555) and U.S. Provisional Patent Application No. 60/299,129, filed June 18, 2001 (Atty. Docket No. MBHB00-555-A). All patents, patent applications (published or unpublished) and other scientific or technical writings referred to herein are hereby incorporated by reference to the extent that they are not contradictory.

BACKGROUND OF THE INVENTION

15

Field of the Invention

The invention is related to the field of blood analysis. In particular, blood analysis using a solid phase coated with marker-specific compounds. More in particular, blood analysis wherein the solid phase comprises paramagnetic particles and the marker-specific compounds comprise antibodies, receptors, ligands, proteins, peptides, cytokines, chemokines, small molecules and the like. Even more in particular, the invention is related to the analysis of activated and unactivated platelets from whole blood as well as the identification of chemical markers associated with the coagulation process.

25

Description of Related Art

Separation of platelets from whole blood is generally accomplished by centrifugation of the blood at 150 x g for 10 min. The platelet-rich plasma (PRP) fraction is then carefully removed from the top layer and the platelets are subsequently isolated for later use such as transfusion and/or diagnostic tests of platelet function. However, platelets are very prone to artifactual activation when they are centrifuged and handled (Metcalfe, P., Williamson, L.M., Reutelingsperger, C.P., Swann, I., Ouwehand, W.H. & Goodall, A.H. (1997) "Activation during preparation of therapeutic platelets affects deterioration during

storage: a comparative flow cytometric study of different production methods." *Br. J. Haematol.*, vol. 98, no. 1, pp. 86-95). Artifactual activation has hampered the platelet diagnostic field for decades as the platelet releases many substances upon activation that could be utilized as new diagnostic markers in the identification of patients with "hyperactive platelets." However, often platelet specific markers, such as beta-thromboglobulin (Mundal, H.H., Hjemdahl, P., Urdal, P., Kierulf, P., Perneby, C., Bergt, K. & Gjesdal, K. (1998) "Beta-thromboglobulin in urine and plasma: influence of coronary risk factors." *Thromb. Res.*, vol. 90, no. 5, pp. 229-237) and thromboxane B-2 (Ciabattini, G., Maclouf, J., Catella, F., FitzGerald, G.A. & Patrono, C. (1987) "Radioimmunoassay of 11-dehydrothromboxane B2 in human plasma and urine." *Biochim. Biophys. Acta*, vol. 918, no. 3, pp. 293-297) are increased as a consequence of isolating the platelets for assay by centrifugation and further handling in the assay. Thus, these markers have not been adapted in the identification of platelet activation as initially hoped. It is an object of this invention to provide a novel method for separating platelets from whole blood, without centrifugation, such that subsequent analysis of platelet specific markers can be accomplished without artifactual elevation of such markers due to processing.

Platelet activation and subsequent aggregation are known to play a pivotal role in the acute pathophysiology of thrombus formation, stroke and acute coronary syndromes (ACS). ACS patients experiencing unstable angina and/or non Q-wave myocardial infarction are prone to plaque rupture and thrombus formation which is amenable to a host of pharmacological agents (thrombolytics, GPIIb/IIIa antagonists, anti-coagulants). Thus, it is imperative that platelet activation be assessed rapidly and accurately such that appropriate therapeutic interventions be employed to salvage ischemic myocardium that is at risk of infarction (death). Methods that allow for separation of platelets from whole blood without centrifugation, will be valuable in developing new assays for detecting "hyperactive" platelets that may contribute to disease states.

All platelets, either unactivated, activated or circulating as microparticles, express glycoprotein 1b (GP1b) on their surface (White, J.G., Krumwiede, M.D. &

Escolar, G. (1999) "Glycoprotein 1b is homogeneously distributed on external and internal membranes of resting platelets." *Am. J. Pathol.*, vol. 155, no. 6, pp. 2127-2134). Thus, GP1b represents an appropriate target for identifying all circulating platelets and platelet microparticles. In the present invention, anti-
5 GP1b monoclonal antibodies may be employed to "capture" platelets in which the antibody is coated onto the surface of paramagnetic particles. Mixing whole human blood with GP1b-coated paramagnetic particles, followed by magnetic separation, results in platelet capture and substantial depletion of platelets from a given sample.

10

BRIEF SUMMARY OF THE INVENTION

This invention relates to the analysis of biological samples using one or more coated phase(s). The biological sample may be, but is not limited to,
15 undiluted and/or diluted whole blood, undiluted and/or diluted blood plasma, as well as given fraction(s) of fractionated whole blood. The solid phase may be, but is not limited to, paramagnetic particles. The solid phase may be coated with one or more marker-specific protein tracers including, but not limited to, antibodies, receptors, ligands, proteins, peptides, cytokines, chemokines, small
20 molecules and the like.

This invention further relates to methods for separating platelets without activating the platelets during the separation process to provide a platelet sample containing activated and unactivated platelets along with platelet-derived
25 microparticles wherein the activated platelets are activated by physiological processes *in vivo* and not by the separation process.

The invention includes separated platelet compositions that are substantially unactivated by the separation process and comprising both
30 physiologically activated platelets and unactivated platelets along with platelet-derived microparticles. The invention includes assays of the separated platelet samples based on the markers from the physiologically activated platelets.

A preferred method for separating platelets and platelet derived microparticles involves attaching an antibody or protein that specifically binds to platelets and platelet derived microparticles onto paramagnetic particles and contacting a diluted or undiluted whole blood sample or fraction thereof with the antibody-coated paramagnetic particles, magnetically separating the paramagnetic particles and attached platelets and removing the remaining substantially platelet free supernatant from the paramagnetic particles-platelet complexes. This separation provides a platelet composition where platelets are not activated by the separation process and the only activated platelets in the composition are those that have been physiologically activated *in vivo*. Thus, analysis of cellular markers in this composition is a more accurate measurement of the *in vivo* physiological activated platelets.

This invention also relates to methods for separating platelets and platelet-derived microparticles without activating the platelets during the separation process to provide a sample substantially free of platelets and platelet-derived microparticles.

The invention includes assays of the substantially platelet-free samples based on markers which could be significantly influenced by the presence of physiologically-activated platelets within the sample.

A preferred method for obtaining a substantially platelet free sample involves separating platelets by attaching an antibody that specifically binds to platelets onto paramagnetic particles and contacting a whole blood sample or diluted whole blood sample with the antibody-coated paramagnetic particles, magnetically separating the paramagnetic particles and attached platelets and separating the remaining platelet-free supernatant from the paramagnetic particle/platelet complexes. This separation process provides a sample composition substantially free of platelets and physiologically unaltered by the separation process. Thus, analysis of soluble markers in this composition is a more accurate measurement of the *in vivo* physiological state.

It is also contemplated that this invention would be applicable to the identification of other diseases and conditions. By judicious selection of the material coating the paramagnetic particles (e.g., including, but not limited to, antibodies (polyclonal and/or monoclonal), ligands, receptors, proteins, peptides, cytokines, chemokines, small molecules and the like), one would be able to analyze undiluted and/or diluted whole blood, and/or any fraction thereof, for the presence or absence of markers of these other diseases and conditions. These other diseases and conditions, and representative markers, include, but are not limited to, the following:

Atherosclerosis/Inflammation/Arthritis:

Proteins: VCAM-1, ICAM-1, P- and E-selectin, P-selectin glycoprotein ligand-1 (PSGL-1), PECAM-1, C-reactive protein (hCRP), ox-LDL, HDL, LDL, Apolipoprotein A1 (Apo A-1), total cholesterol, LP(a), CD15, CD40, Interleukin-6 (IL-6), Interleukin-1 receptor antagonist (IL-1ra), Tumor Necrosis Factor (TNF), Tissue Factor (TF), Tissue Factor Pathway Inhibitor (TFPI), Complement C3a and C5a, C3 and C5 Convertase, Factor D, Kallikrein, Plasmin, C1-Inhibitor, soluble CR1, etc.

Chemokines: Monocyte Chemoattractant Protein-1 (MCP-1), MCP-4, Regulation on Activation Normal T-cell Expressed and Secreted (RANTES), Interleukin-8 (IL-8), Stromal Cell Derived Factor-1 (SDF-1), etc.

Diabetes:

Proteins: P- and E-selectin, GPIIb/IIIa, lysosomal GP53, thrombospondin, Glucose associated Hemoglobin (GHb), glycohemoglobin A(1c), gamma globulin, Insulin-like growth factor (IGF-1), Insulin-like Growth Factor Binding Proteins 1-6, Pregnancy Associated Plasma Protein A (PAPP-A), Receptor for Advanced Glycation End-Products (RAGE), etc.

Transplant Rejection:

Proteins: hCRP, E- and L-selectin, CD18, CD11a, CD11b, Interleukins (IL-1, IL-2, IL-6, IL-15), Complement C3a and C5a, C3 and C5 Convertase, Factor D, soluble CR1, Interferon gamma (IFN-gamma), Chemokines (CCR2, CCR3, CCR5), TNF alpha, IgG, etc.

Alzheimer's Disease/Vascular Dementia:

Proteins: Beta-amyloid peptide, Protein tau, Beta-amyloid precursor protein, Abeta1-40, Abeta1-42, presenillins, ApoE4, C3a, C5a, C3 and C5 Convertase, ERp57, etc.

Additional Thrombosis/Stroke/Peripheral Arterial Disease Targets (other than those already mentioned (i.e., GPIIb/IIIa, P-selectin, D-dimer, PTF1.2)):

Proteins: Factor (FXa), Factor VIIa (FVIIa), Factor IXa (FIXa), Factor V(FV), Thrombin (FIIa), Factor XIIIa (FXIIIa), Tissue Factor (TF), Tissue Factor Pathway Inhibitor (TFPI), CD45, activated Partial Thromboplastin Time (aPTT), Prothrombin Time (PT/INR), Thrombomodulin, Thrombospondin, Thromboxane (TxA-2), Plasminogen Activator Inhibitor-1 (PAI-1), Thrombin Activatable Fibrinolysis Inhibitor (TAFI), Tyrosine Kinase SYK, Angiotensin IV, P2T, von Willebrand's Factor (vWF), Fibrinopeptide A (FPA), Fibrinopeptide B (FPB), Fibrin Degradation Products (FDP's), Thrombin-Antithrombin Complex (TAT), Pentraxin (PTX3), t-PA (tissue plasminogen activator), u-PA (urokinase plasminogen activator), Plasminogen, Plasmin, Factor XIII, alpha-2-plasmin inhibitor, alpha-1 anti-trypsin inhibitor, etc.

Hemophilia and Related Disorders:

Proteins: Factor VIII, Factor IX, etc.

Others:

Proteins: peptides and proteins generated during platelet activation and coagulation, neurological and neuro-associated peptides and proteins, antibodies to heparin in heparin associated thrombocytopenia, etc.

5

BRIEF DESCRIPTION OF THE SEVERAL VIEWS OF THE DRAWINGS

Figure 1 shows a graph illustrating correlation of platelet counts of whole blood diluted in Cellpack diluent and in phosphate buffered saline (PBS), pH 7.4
10 supplemented with 5% bovine serum albumin (BSA) and 10% CTAD.

Figure 2 shows a graph illustrating a linear response of platelet counts in diluted whole blood.

15 Figure 3 shows the percent platelet capture versus capture time in diluted whole blood. The capture data for three concentrations of coated particles is displayed.

Figure 4 shows a correlation of fluorescence intensity versus soluble P-selectin (at various concentrations in diluted plasma) captured using anti-P-selectin (anti-CD62P)-coated paramagnetic microparticles.
20

Figure 5 shows a correlation of fluorescence intensity versus membrane P-Selectin for different sample preparation conditions using a 1-step assay
25 format.

Figure 6 shows a correlation of fluorescence intensity versus membrane P-Selectin for different sample preparation conditions using a 2-step assay
format.
30

Figure 7 shows a correlation of fluorescence intensity versus membrane GPIIb/IIIa for different sample preparation conditions (*i.e.*, sample volume).

Figure 8 shows a correlation of fluorescence intensity versus membrane P-Selectin for different sample preparation conditions (*i.e.*, sample volume).

DETAILED DESCRIPTION OF THE INVENTION

5

The term "label" means a group or compound attached to an antibody or an analyte or an analyte analogue that renders the reaction between the antibody or analyte or analyte analogue detectable. Representative examples of labels include enzymes, radioactive elements, fluorophores, and chemicals that produce light. A label is any substance, either alone or in conjunction with other substances, that can be attached to an appropriate molecule and that is capable of producing a signal that is detectable by visual or instrument means. Various labels include catalysts, enzymes, liposomes, and other vesicles containing signal producing substances such as chromogens, catalysts, fluorescent compounds, chemiluminescent compounds, enzymes, radioactive elements and the like. In this invention, the preferred label is fluorescent. The term "tracer" is synonymous with the term "label".

The term "solid phase" means a plurality of microparticles having specific binding members chemically or physically bound thereto. Other solid phases that are known to those skilled in the art include the walls of wells or reaction trays, tubes, polymeric beads, nitrocellulose strips, membranes, chromatographic columns and the like. A preferred solid phase comprises microparticles made of polystyrene containing a layer of iron oxide rendering them paramagnetic. The preferred method of separating the particles from the test sample involves capture of the particles by means of a magnetic field. In this invention the preferred method, the solid phase consists of paramagnetic microparticles having specific binding members chemically or physically bound thereto.

The term "sample", "biological sample", and the like mean a material suspected of containing an analyte. The sample can be used directly as obtained from the source or following a pretreatment to modify the character of the sample. The sample can be derived from whole blood. The sample can be

treated prior to use, such as preparing plasma from whole blood, diluting viscous fluids, and the like. Methods of treatment can include fractionation, filtration, extraction, concentration, the addition of reagents and the like.

5 The following non-limiting examples will further explain the invention

Example 1

Platelet Capture

10 This example illustrates the capture and removal of platelets from whole blood or diluted whole blood.

Materials and Methods

15 Whole blood was obtained from healthy volunteers and collected into centrifuge tubes containing citrate theophylline adenosine dipyridamole (CTAD), pH 5.4 anticoagulant at a 9:1 ratio. A volume of whole blood was diluted to 1.85% in phosphate buffered saline (PBS), pH 7.4, supplemented with 5% bovine serum albumin (BSA) and 10% CTAD.

20

Platelet cell counts were examined by means of an automatic analyzer (Sysmex Microcellcounter F-800 Hematology Analyzer, Sysmex Corp. of America, Long Grove, IL; see also Fujimoto, K. (1999) "Principles of Measurement in Hematology Analyzers Manufactured by Sysmex Corporation" *Sysmex Journal International*, vol. 9, no. 1, pp. 31-44). Each sample condition was pipetted into a Sysmex Disposable Sample Beaker (DB-1) and further diluted with 10.0 mL of Cellpack Whole Blood Diluent (CPK-310A) and assayed in duplicate.

25

30 Platelet counts from whole blood were compared as a function of dilution with either Cellpack (buffer supplied with the Sysmex Microcellcounter F-800 Hematology Analyzer system) or with phosphate-buffered saline supplemented with 5% bovine serum albumin (PBS/BSA) diluent. The data shown in Figure 1

demonstrates that platelet counts decreased as a function of whole blood dilution and platelet counts obtained from (PBS/BSA) diluent correlated well with the Sysmex Hematology Analyzer Cellpack diluent.

5 In order to demonstrate measurements were performed within the linear detection portion of the Sysmex Microcellcounter F-800 Hematology Analyzer, the following experiments were conducted: 20 uL of whole blood was pipetted into 980 uL of PBS containing 5% BSA, 10% CTAD and then serially dilute as appropriate. 100 uL of prediluted sample were deposited into a DB-1 Sample
10 Beaker and further diluted with 10.0 mL of Cellpack Whole Blood Diluent for Platelet cell count measurements. The data shown in Figure 2 demonstrates that platelet counts in 1.85%-diluted were performed within the linear portion of the Sysmex platelet detection curve. Plotted is a dilution curve from 0-2% with a corresponding high correlation coefficient (0.9973).

15

3.7 micron (μm) carboxyl-modified paramagnetic microparticles were coated with murine monoclonal antibody (mAb) which specifically recognizes the platelet membrane surface component glycoprotein GP1b, (CD42b, Biodesign, N42409M, lot 6G1996) and were re-suspended in phosphate buffered saline
20 (PBS)- supplemented with 1% bovine serum albumin (BSA) buffer at 7.4% w/v. 50 μL of the paramagnetic microparticle preparation was added to a reaction tube, separated magnetically and supernatant removed. The particles were resuspended with 100 μL of 1.85% whole blood dilution and incubated at room temperature for 30 minutes. The paramagnetic particles were separated
25 magnetically and supernatant removed for platelet cell count analysis.

Table 1 represents the results of this experiment. Whole human blood was diluted to 1.85% and separation of platelets performed with murine monoclonal antibody anti-GP1b coated paramagnetic microparticles. The
30 samples were incubated for 30 minutes at room temperature and then the microparticles were separated from the solution magnetically and the supernatant removed. The platelet counts from each reaction supernatant was

pipetted into a disposable sample beaker (DB-1) and further diluted with 10 ml of Cellpack diluent and platelet counts determined in the Hematology Analyzer.

Incubation of whole blood with anti-GP1b murine monoclonal antibody
5 coated paramagnetic particles resulted in approximately 90% capture of platelets
in the sample. These results demonstrate that paramagnetic microparticles,
coated with a murine monoclonal antibody directed against platelet GP1b,
followed by magnetic separation, are capable of capturing platelets from human
10 whole blood with no other traditional manipulations such as centrifugation for
subsequent use of the platelet fraction in determining the state of platelet
activation in patients.

Following the above methods, monoclonal antibodies to antigens listed
below and magnetic particles coated with such antibodies are used to separate
15 platelets without activation:

CD42a (GPIX) - part of GP1b-IX-V complex;

CD42b (GP1b-alpha) - 145 kD platelet binding site for vWF and thrombin;

20

CD42c (GP1b-beta) - 25 kD disulfide bonded to alpha subunit;

CD42d (GPV);

25 CD41 (GPIIb, also known as alpha IIB integrin);

CD61 (GPIIIa) - beta 3 subunit of GPIIb/IIIa complex (alpha 2b, beta 3);

30 CD41/CD61 (GPIIb/IIIa complex) - receptor for fibrinogen, fibronectin, von
Willebrand factor, and other adhesion proteins containing the Arg-Gly-Asp motif

CD36 (GPIV) - platelets/monocytes;

CD49b (VLA-2) - platelets/monocytes;

CD51 (alpha V, beta 3) - vitronectin receptor;

5 CD62p (P-selectin) - platelets; and

CD107a (LAMP-2) - lysosomal protein translocated to cell surface after activation

10 CD41a (GPIIb/IIIa) - intact IIb/IIIa complex; fibrinogen, von Willebrand factor, fibronectin and vitronectin receptor.

Table 1. Platelet Capture in Whole Human Blood

Test 1

Results:		<u>Avg.*</u>	<u>Avg. - bkgd.*</u>
No Particle Control	1334, 1302	1318	1309
Positive Control 1**	89, 104	95	86
Positive Control 2**	87, 101		
Sysmex System Buffer	10, 8	9	
Percent Capture	93.4% Capture		

*(cell count x 10^3 cells / uL)

** positive control (anti-GP1b paramagnetic particles)

Test 2

Results:		<u>Avg.*</u>	<u>Avg. - bkgd.*</u>
No Particle Control	1121, 1086,	1111.3	1093.3
No Particle Control	1102, 1136		
Positive Control 1**	163, 159	139.8	121.8
Positive Control 2**	116, 121		
Sysmex System Buffer	17, 19	18.0	
Percent Capture	88.9% Capture		

*(cell count x 10^3 cells / uL)

** positive control (anti-GP1b paramagnetic particles)

Test 3

	<u>Results</u>		<u>Avg.*</u>	<u>Avg. - bkgd.*</u>
5	No Particle Control	1142,1173	1144.5	1126.5
	No Particle Control	1130,1133		
	Positive Control 1**	30, 72	63.0	45.0
10	Positive Control 2**	73, 77		
	Sysmex System Buffer	17, 19	18.0	
	Percent Capture	96.0% Capture		
15	*(cell count x 10 ³ cells / uL)			
	** positive control (anti-GP1b paramagnetic particles)			

20

Summary of Tests 1, 2 & 3
(n=3 human subjects)

		<u>% Capture</u>
25	Test 1	93.4
	Test 2	88.9
	Test 3	96.0
	Mean \pm SEM	92.8 \pm 3.6 %
30		

Example 2

Platelet Capture in Diluted Whole Blood

This example illustrates the dose response of time and solid phase
5 percentage (surface area) to percent platelet capture from whole blood samples.

A multivariant designed experiment was used to examine platelet capture
in diluted whole blood. Experimental parameters were: incubation time was
about 3 - 7 minutes; antibody coating concentration was about 6 - 50 μg ; particle
10 concentration was about 4 - 10 % (w/v). The assay was carried out as follows:
whole blood was diluted to 2% with phosphate buffered saline (PBS) at pH 7.2
supplemented with 1% bovine serum albumin (BSA). 100 μL of paramagnetic
anti-GP1b-coated particles were added to 100 μL of diluted whole blood at room
temperature. Following an incubation period, the particles were separated
15 magnetically and the supernatant was removed for testing in a Sysmex F-800
Hematology Analyzer Microcellcounter as previously described in example 1.
The platelet concentration per μL was determined for each condition and
compared to the original 2% stock solution. Particle coating concentration,
particle % (w/v) and incubation times were varied. The results of these
20 experiments are displayed in Figure 3. The data demonstrates that effective
platelet removal from whole blood is effected significantly by the percentage of
solid phase utilized. Through the addition of increased higher percent solid
phase concentrations, in turn ever increasing percentage of platelets are
removed from the sample. The results indicate that for any given whole blood
25 sample or whole blood sample dilution effective platelet capture is achievable if
sufficient time, antibody coating concentration and percentage of anti-platelet
coated solid phase is utilized.

Example 3

30

This example further illustrates the capture and removal of platelets from
whole blood or diluted whole blood.

Whole blood was obtained from healthy volunteers and collected into centrifuge tubes containing D-phe-pro-arg-chloromethylketone (PPACK, dihydrochloride) anticoagulant. A volume of whole blood was diluted to 2.0 % in phosphate buffered saline (PBS), pH 7.4, supplemented with 1% bovine serum albumin (BSA)

100 uL of Murine monoclonal antibody anti-GP1b coated paramagnetic microparticles were added to 100 μ L of a 2.0% whole blood dilution and incubated at room temperature for 5 minutes. The paramagnetic particles were separated magnetically and the supernatant removed for platelet cell count analysis as previously described in Example 1 with a Sysmex Microcellcounter F-800 Hematology Analyzer.

Table 2

			Avg.*	Avg. - bkgd.*
15	No Particle Control	601, 605	601	599
	No Particle Control	594, 605		
	Positive Control 1**	67, 68	59	57
20	Positive Control 2**	52, 55		
	Positive Control 3**	53, 58		
	Sysmex System Buffer	2, 4, 3, 1, 1	2	
25	Percent Capture	90.5% Capture		

*(cell count x 10^3 cells / uL)

** positive control (anti-GP1b paramagnetic particles)

30

Example 4

Soluble P-selectin Assay in Plasma

5 This example illustrates the feasibility of assaying substantially platelet free samples for markers which can be influenced by the presence of physiologically activated platelets within the sample. In this example soluble P-selectin concentrations are measured in substantially platelet free samples to allow for discrimination from membrane P-selectin.

10 Recombinant human P-selectin was added to 2.0% platelet-free human plasma sample at the following concentrations: 400, 200, 100, 50, 25, 10 and 1.0 ng/mL. 50 μ L of 0.125% (w/v) paramagnetic microparticles (1.5 micron (μ m)) coated with an anti-P-selectin murine monoclonal antibody (Mab) (anti-CD62P)
15 were added to 100 μ L of each of the aforementioned sample concentrations of recombinant P-selectin in plasma at room temperature and incubated for 10 minutes. The microparticles were separated magnetically and the supernatant was removed. The microparticles were washed by re-suspension with 200 μ L of PBS supplemented with 0.5% (BGG) diluent. The microparticles were separated
20 magnetically and the supernatant was removed. The washed microparticles were re-suspended in 50 μ L of a 20 μ g/mL solution of FITC-labeled anti-P-selectin (anti-CD62P) rabbit polyclonal antibody and incubated at room temperature for 10 minutes. The microparticles were separated magnetically and the supernatant removed. The microparticles were washed by re-suspension
25 with 200 μ L of (PBS) supplemented with 0.5% BGG diluent. The wash sequence was repeated and the particles were re-suspended in 0.1% SDS at pH 11. The microparticles were separated magnetically and the supernatants were transferred to a 96 well microtiter-plate and examined for fluorescence intensity. The graph in Figure 4 illustrates how the incubation of samples with
30 paramagnetic microparticles coated with anti-P-selectin murine monoclonal antibody followed by the addition of FITC labeled anti-P-selectin rabbit polyclonal antibody enables the detection of and quantitation of soluble P-selectin.

Example 5
1-Step Assay

Whole blood was obtained from a healthy volunteer and collected into a blood collection tube containing D-phe-pro-arg-chloromethylketone (PPACK, Dihydrochloride). An aliquot of whole blood was stimulated with 10 μ M adenosine diphosphate (ADP) for two minutes and then fixed for two minutes with neutral buffered formalin (1% final concentration). Another aliquot was not treated with ADP, but only fixed with 1% neutral buffered formalin for two minutes

To 100 μ L of each of the above three sample preparations the following reagents were added: 50 μ L of FITC-labeled anti-P-selectin (anti-CD62P) rabbit polyclonal antibody (at 20 μ g/mL) and 50 μ L of 0.5% (w/v) paramagnetic microparticles (3.7 micron (μ m)) coated with anti-glycoprotein GP1b (CD42b) murine monoclonal antibody (Mab). The mixtures were incubated at room temperature for either 10 or 15 minutes. The microparticles were separated magnetically and the supernatant was removed. The microparticles were washed by re-suspension in 200 μ L of PBS supplemented with 1% bovine serum albumin (BSA) followed by magnetic separation of the microparticles and removal of the supernatant. This sequence was repeated and the microparticles were re-suspended in 0.1% SDS at pH 11. The microparticles were separated magnetically and the supernatants were transferred to a 96 well plate and examined for fluorescence intensity.

The data in table 3 as well as the graph in Figure 5 illustrates how the incubation of whole blood with both paramagnetic microparticles coated with an anti-GP1b murine monoclonal antibody and a labeled anti-P-selectin rabbit polyclonal antibody enables the detection and discrimination of platelets that are expressing membrane-bound P-selectin from those that are not expressing membrane-bound P-selectin.

Table 3

FITC Signal Detected vs. Sample Preparation Conditions
(at representative incubation times)

<u>Sample Preparation Conditions</u>	<u>10 min.</u>	<u>15 min.</u>
ADP agonist treated whole blood	1235	1526
Non agonist treated whole blood	304	282

Example 6
2-Step Assay

5

Whole blood was obtained from a healthy volunteer and collected into a sample collection tube containing D-Phe-Pro-Arg-chloromethylketone (PPACK, dihydrochloride). An aliquot of whole blood was stimulated with 10 μ M adenosine diphosphate (ADP) and 1.0 μ M Epinephrine for two minutes and then treated for two minutes with neutral buffered formalin (1% final concentration). An additional aliquot of whole blood was stimulated with 10 μ M ADP for two minutes and then treated for two minutes with neutral buffered formalin. An additional aliquot of whole blood was treated for two minutes with 1% neutral buffered formalin for two minutes without ADP stimulation.

15

50 μ L of 0.5% (w/v) paramagnetic microparticles (3.7 micron (μ m)), which were coated with anti-glycoprotein GP1b murine monoclonal antibodies (Mab), were added to 100 μ L of each of the above sample preparations. The mixtures were incubated at room temperature for 10 minutes. The microparticles were separated magnetically and the supernatant was removed. The microparticles were washed by re-suspension in 200 μ L of PBS supplemented with 1% BSA. The microparticles were again separated magnetically and the supernatant removed. The microparticles were re-suspended in 50 μ L FITC-labeled anti-P-selectin (anti-CD62P) rabbit polyclonal antibody at 20 μ g/mL. The assays were incubated at room temperature for an additional 10 minutes. The microparticles were separated magnetically and the supernatant removed. The microparticles

20

25

were washed by re-suspension in 200 μ L PBS supplemented with 1% BSA. The microparticles were again separated magnetically and the supernatant removed. The wash sequence was repeated and the microparticles re-suspended in 0.1% SDS at pH 11. The microparticles were separated magnetically and the supernatants transferred to a 96 well plate and examined for fluorescence intensity.

The data in table 4 as well as the graph in Figure 6 illustrate how the incubation of whole blood with paramagnetic microparticles coated with anti-GP1b murine monoclonal antibodies followed by incubation with labeled anti-P-selectin rabbit polyclonal antibodies enables the detection, isolation and discrimination of platelets that express membrane P-selectin from those platelets that do not express membrane P-selectin.

Table 4

FITC Signal Detected vs. Sample Preparation Conditions

<u>Sample Preparation Conditions</u>	<u>FITC signal detected</u>
ADP & Epinephrine agonist treated whole blood	19550
ADP agonist treated whole blood	17349
Non agonist treated whole blood	2986

Example 7

GPIIb/IIIa Assay

3.7 micron (μm) carboxyl-modified paramagnetic microparticles were coated with Human Fibrinogen which interacts and binds with the complexed form of the platelet membrane surface component glycoproteins GP1Ib/GPIIIa (CD41/CD61). The coated particles were re-suspended in phosphate buffered saline (PBS) buffer at 0.5% w/v.

Whole blood was obtained from a healthy volunteer and collected into a 3.0 mL blood collection tube containing 0.3 μg D-phe-pro-arg-chloromethylketone (PPACK, Dihydrochloride) and 2.1 mL Hanks' Balanced Salt supplemented with 10 mM HEPES, pH 7.4. An aliquot of this whole blood preparation was stimulated with 10 μM adenosine 5' -diphosphate (ADP) for two minutes prior to being tested in the assay. Another aliquot was stimulated with 2 μM adenosine 5' -diphosphate (ADP) for two minutes prior to being tested in the assay, and a final aliquot was not treated with ADP.

15

To a set of 5 μL , 10 μL , and 20 μL of each of the above three sample preparations 20 μL of FITC-labeled anti-IIb/IIIa (anti-CD41a) mouse monoclonal antibody (at 20 $\mu\text{g/mL}$) was added. The mixture was incubated for 10 minutes at room temperature at which point 50 μL of 0.5% (w/v) paramagnetic microparticles (3.7 micron (μm)) coated with Fibrinogen was added. The mixtures were incubated at room temperature for an additional 10 minutes. The microparticles were separated magnetically and the supernatant was removed. The microparticles were washed by re-suspension in 200 μL of PBS supplemented with 1% bovine serum albumin (BSA), followed by magnetic separation of the microparticles and removal of the supernatant. This sequence was repeated three times and the microparticles were re-suspended in 0.1% SDS at pH 11. The microparticles were separated magnetically and the supernatants were transferred to a 96 well plate and examined for fluorescence intensity.

30

The data in table 5, as well as the graph in Figure 7, illustrate how the incubation of whole blood with both paramagnetic microparticles coated with fibrinogen and a labeled anti-GPIIb/IIIa monoclonal antibody enables the detection and discrimination of whole blood samples containing platelets that are

expressing membrane-bound GPIIa/IIIb from those that are not expressing membrane-bound GPIIa/IIIb.

Table 5

FITC Signal Detected vs. Sample Volume

	<u>5 μL</u>	<u>10 μL</u>	<u>20 μL</u>
10 μ M ADP agonist treated whole blood aliquot	525	616	1031
2 μ M ADP agonist treated whole blood aliquot	524	579	896
Non agonist treated whole blood aliquot	171	0	98

5

Example 8

Membrane P-Selectin Assay

10 Whole blood was obtained from a healthy volunteer and collected into a 3.0 mL blood collection tube containing 0.3 μ g D-phe-pro-arg-chloromethylketone (PPACK, Dihydrochloride) and 2.1 mL Hanks' Balanced Salt supplemented with 10 mM HEPES, pH 7.4. An aliquot of this whole blood preparation was stimulated with 10 μ M adenosine 5' -diphosphate (ADP) for two minutes prior to
15 being tested in the assay. Another aliquot was stimulated with 2 μ M adenosine 5' -diphosphate (ADP) for two minutes prior to being tested in the assay, and a final aliquot was not treated with ADP.

To a set of 10 μ L and 20 μ L of each of the above three sample
20 preparations 20 μ L of FITC-labeled anti P-Selectin (anti-CD62P) rabbit polyclonal antibody (at 20 μ g/mL) was added. The mixture was incubated for 10 minutes at room temperate at which point 50 μ L of 0.5% (w/v) paramagnetic microparticles (3.7 micron (μ m)) coated with anti-glycoprotein GP1b (CD42b) murine monoclonal antibody (Mab) was added. The mixtures were incubated at room
25 temperature for an additional 10 minutes. The microparticles were separated magnetically and the supernatant was removed. The microparticles were

washed by re-suspension in 200 μ L of PBS supplemented with 1% bovine serum albumin (BSA) followed by magnetic separation of the microparticles and removal of the supernatant. This sequence was repeated three times and the microparticles were re-suspended in 0.1% SDS at pH 11. The microparticles were separated magnetically, and the supernatants were transferred to a 96 well plate and examined for fluorescence intensity.

The data in table 6, as well as the graph in Figure 8, illustrate how the incubation of whole blood with both paramagnetic microparticles coated with an anti-GP1b murine monoclonal antibody and a labeled anti P-selectin rabbit polyclonal antibody, enables the detection and discrimination of whole blood samples containing platelets that are expressing membrane-bound P-Selectin from those that are not expressing membrane-bound P-Selectin.

Table 3

FITC Signal Detected vs. Sample Volume

	<u>10 μL</u>	<u>20 μL</u>
10 μ M ADP agonist treated whole blood aliquot	1749	2747
2 μ M ADP agonist treated whole blood aliquot	1544	2390
Non agonist treated whole blood aliquot	170	219

The above set-out examples serve to illustrate the invention and are not intended to limit it in spirit or scope. Those of skill in the art will appreciate that numerous deviations from the examples described above can be made or performed without exceeding the scope of the invention. All patents and patent applications (published or unpublished) as well as other scientific and technical literature referred to herein are expressly incorporated herein by reference to the extent that they are not contradictory.

CLAIMS

What is claimed is:

- 5 1. A method for analyzing a biological sample, wherein the method comprises the steps of:
- (a) combining the biological sample, a coated solid phase and one or more marker-specific, labeled compounds to form a reaction mixture, and
- 10 (b) analyzing the mixture for the presence and/or concentration of the marker.
2. The method of claim 1, wherein the solid phase is separated from the
- 15 mixture prior to analysis.
3. The method of claim 1, wherein the biological sample comprises undiluted or diluted whole blood.
- 20 4. The method of claim 1, wherein the biological sample comprises undiluted or diluted blood plasma.
5. The method of claim 1, wherein the biological sample comprises a undiluted or diluted fraction of fractionated whole blood.
- 25 6. The method of claim 1, wherein the solid phase comprises paramagnetic particles.
7. The method of claim 6, wherein the solid phase is separated from the
- 30 biological sample with a magnet prior to analysis.
8. The method of claim 6, wherein the particles are coated with antibodies or fragments of antibodies with intact complementary determining regions.

9. The method of claim 8, wherein the antibodies or fragments of antibodies are monoclonal or polyclonal.

10. The method of claim 1, wherein the solid phase is coated with one or more substances selected from the group consisting of polyclonal antibodies, monoclonal antibodies, receptors, ligands, proteins, peptides, cytokines, chemokines, small molecules and fragments of any of the preceding.

11. The method of claim 1, wherein the marker-specific, labeled compounds comprise antibodies.

12. The method of claim 1, wherein the label comprises a fluorophore.

13. The method of claim 1, wherein the marker is selected from the group consisting of: CD42c (GP1b-beta)-25 kD disulfide bonded to alpha subunit; CD42d (GPV); CD41 (GPIIb also known as alpha IIB integrin); CD61 (GPIIIa)-beta 3 subunit of GPIIb/IIIa complex (alpha 2b, beta 3); CD41/CD61 (GPIIb/IIIa complex) - receptor for fibrinogen, fibronectin, von Willebrand factor, and other adhesion proteins containing the Arg-Gly-Asp motif; CD36 (GPIV)-platelets/monocytes; CD49b (VLA-2)-platelets/monocytes; CD51 (alpha V, beta 3)-vitronectin receptor; CD62p (P-selectin)-platelets; CD107a (LAMP-2)-lysosomal protein translocated to cell surface after activation and CD41a (GPIIb/IIIa) - intact IIB/IIIa complex; fibrinogen, von Willebrand factor, fibronectin and vitronectin receptor.

14. The method of claim 1, wherein said solid phase is coated with a substance specific for a marker selected from the group consisting of: CD42c (GP1b-beta)-25 kD disulfide bonded to alpha subunit; CD42d (GPV); CD41 (GPIIb also known as alpha IIB integrin); CD61 (GPIIIa)-beta 3 subunit of GPIIb/IIIa complex (alpha 2b, beta 3); CD41/CD61 (GPIIb/IIIa complex) - receptor for fibrinogen, fibronectin, von Willebrand factor, and other adhesion proteins containing the Arg-Gly-Asp motif; CD36 (GPIV) platelets/monocytes;

CD49b (VLA-2)-platelets/monocytes; CD51 (alpha V, beta 3)-vitronectin receptor; CD62p (P-selectin)-platelets; CD107a (LAMP-2)-lysosomal protein translocated to cell surface after activation and CD41a (GPIIb/IIIa) - intact IIb/IIIa complex; fibrinogen, von Willebrand factor, fibronectin and vitronectin receptor.

5

15. The method of claim 1, wherein said marker-specific, labeled compound is specific for a marker selected from the group consisting of: CD42c (GP1b-beta)-25 kD disulfide bonded to alpha subunit; CD42d (GPV); CD41 (GPIIb also known as alpha IIB integrin); CD61 (GPIIIa)-beta 3 subunit of
10 GPIIb/IIIa complex (alpha 2b, beta 3); CD41/CD61 (GPIIb/IIIa complex) - receptor for fibrinogen, fibronectin, von Willebrand factor, and other adhesion proteins containing the Arg-Gly-Asp motif; CD36 (GPIV)-platelets/monocytes; CD49b (VLA-2)-platelets/monocytes; CD51 (alpha V, beta 3)-vitronectin receptor; CD62p (P-selectin)-platelets; CD107a (LAMP-2)-lysosomal protein translocated
15 to cell surface after activation and CD41a (GPIIb/IIIa) - intact IIb/IIIa complex; fibrinogen, von Willebrand factor, fibronectin and vitronectin receptor.

16. The method of claim 1, wherein a buffer is included in the reaction mixture.

20


17. The method of claim 9, wherein the buffer comprises Hepes buffer, Hank's balanced salts and PPACK.

18. The method of claim 2, wherein the biological sample comprises
25 undiluted or diluted whole blood; the solid phase comprises paramagnetic particles; the paramagnetic particles are coated with one or more polyclonal antibodies, one or more monoclonal antibodies, and/or fibrinogen; the marker-specific, labeled compound comprises one or more fluorophore-labeled polyclonal and/or monoclonal antibodies that recognize activated and/or
30 unactivated platelets; the analysis measures the presence, absence and/or concentration of activated and/or unactivated platelets; and the solid phase is separated magnetically from the biological sample.

19. Cells separated from whole blood according to the method of claim 18, wherein the separated cells comprise physiologically activated platelets, unactivated platelets and/or platelet-derived microparticles.

- 5 20. A method for analyzing the biological sample and separated cells of claim 18 in order to assess the relative risk of acute coronary syndrome(s) in the patient or subject from whom the biological sample was obtained.

Figure 1

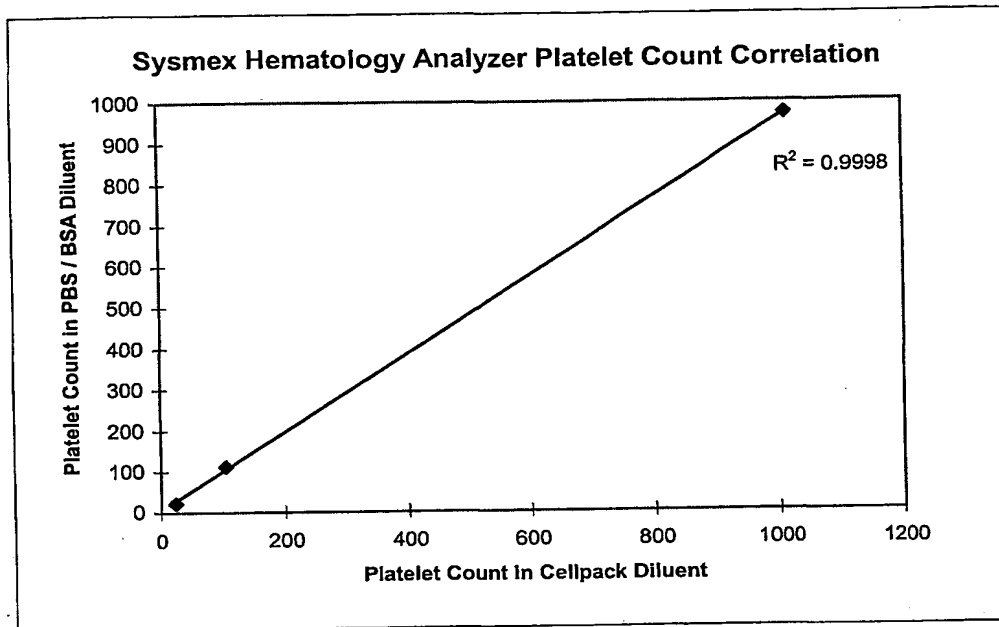

<u>Dilution Factor</u>	<u>Platelet Cell Counts*</u>	
	<u>Cellpack</u>	<u>PBS/BSA</u>
5000	1013	970
50000	105	112
500000	24	22


Note: Platelet Cell Counts $\times 10^3$ cells per microliter

2/8

Figure 2

% Whole Blood Platelet Count

2.0	1421
1.0	781
0.5	379
0.0	7

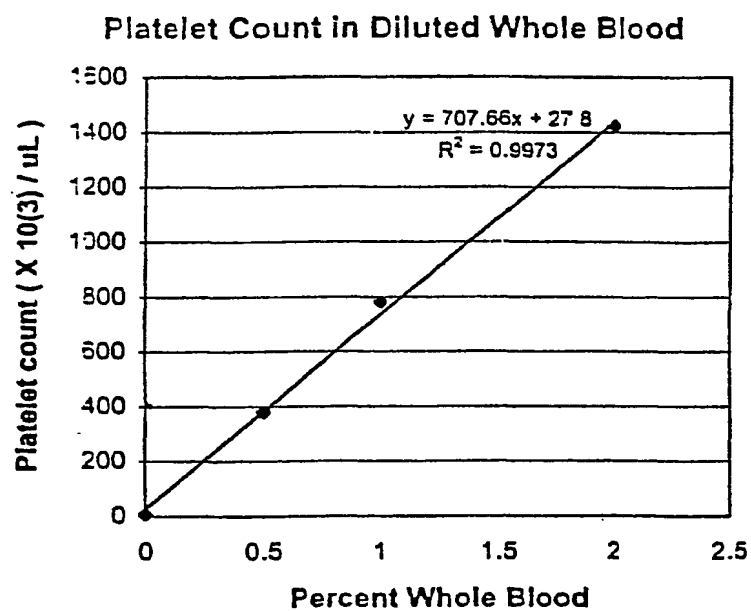

3/8

Figure 3



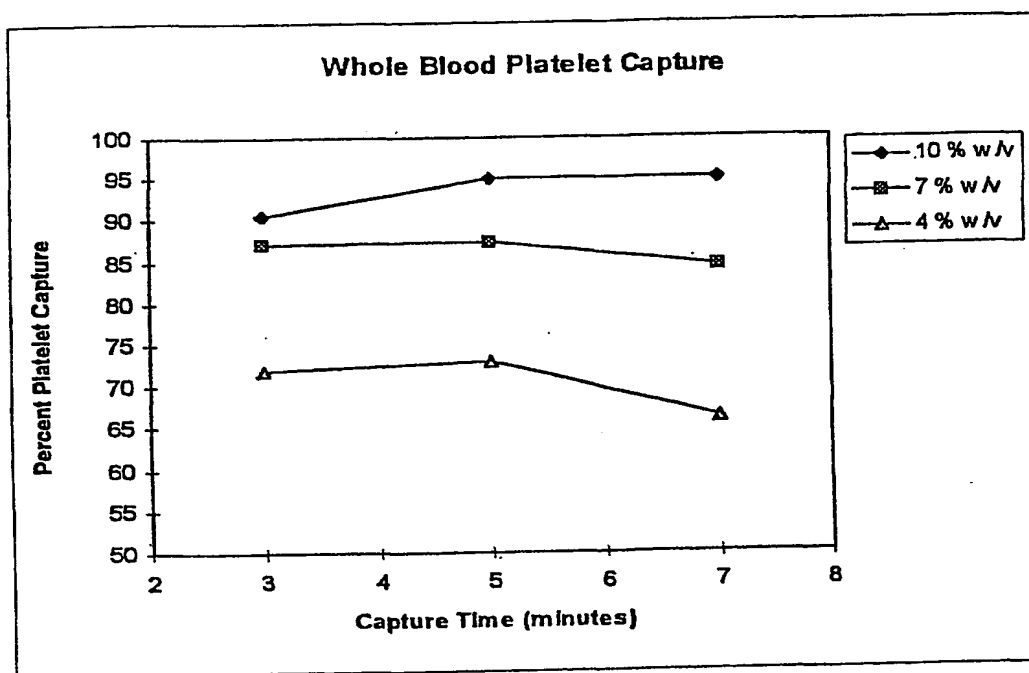


Figure 4

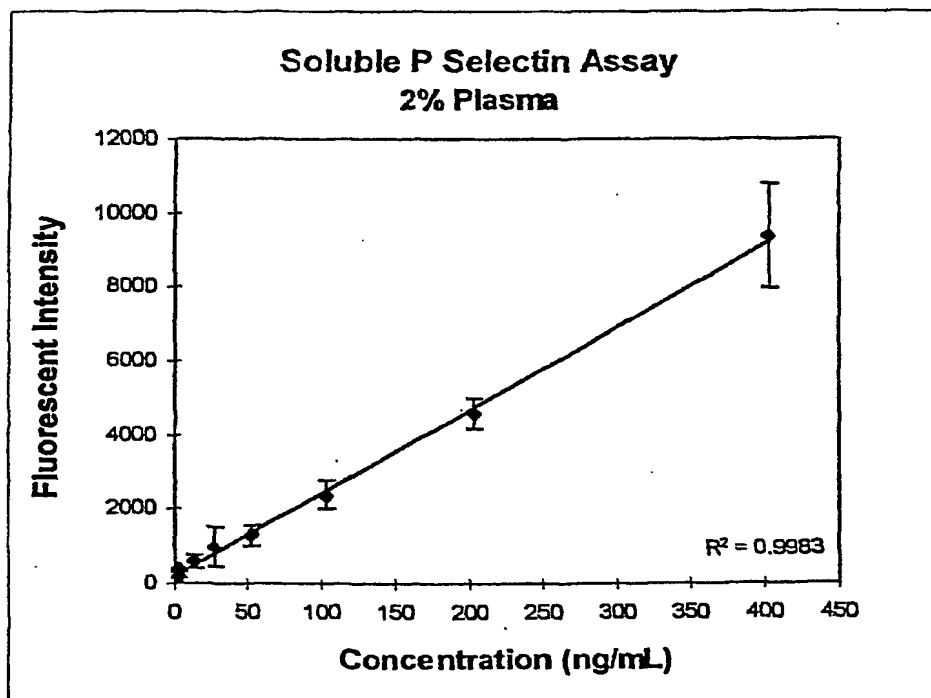


5/8

Figure 5

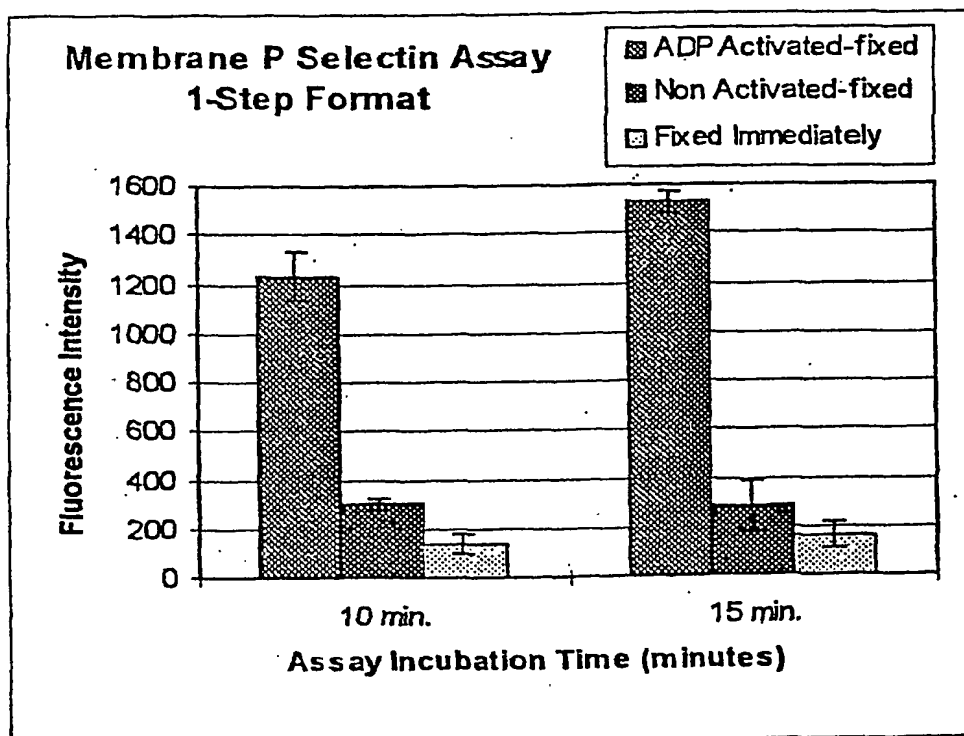


6/8

Figure 6

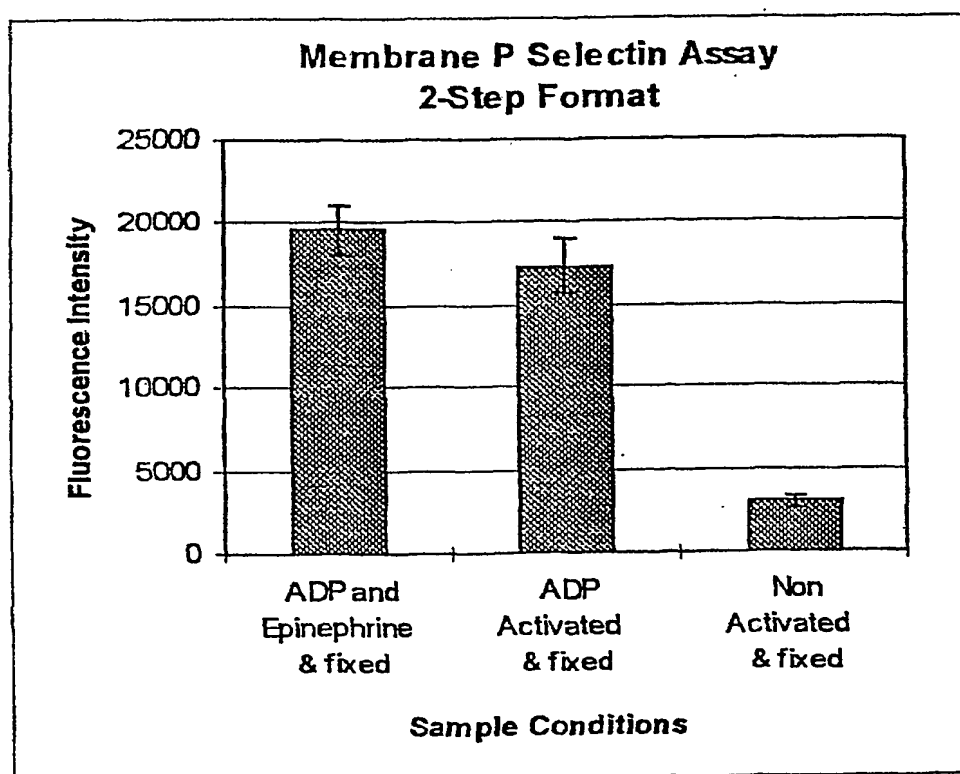


7 / 8

Figure 7

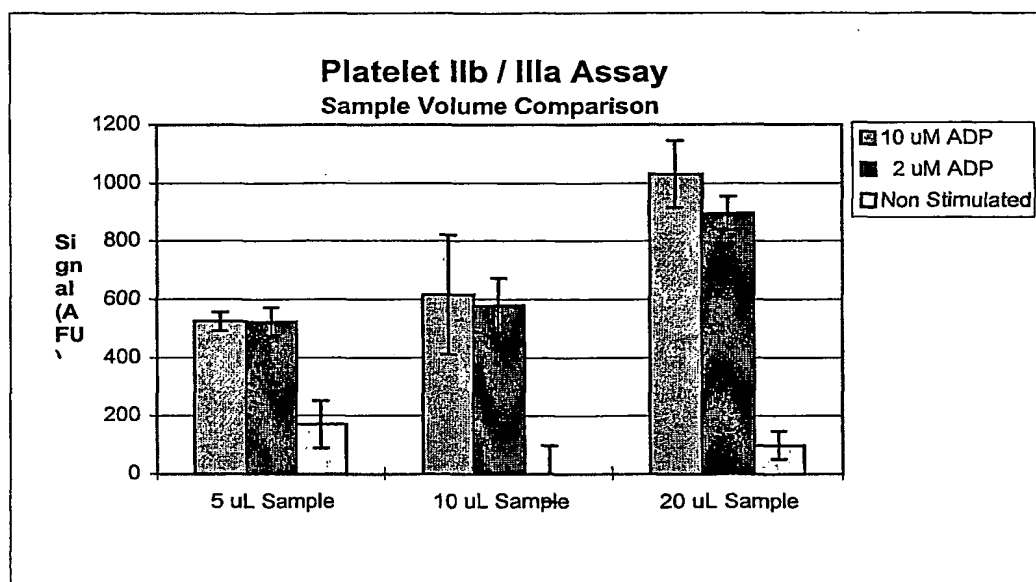

	Fibrinogen Anti-41a	Fibrinogen Anti-41a	Fibrinogen Anti-41a
	5 uL Sample	10 uL	20 uL
10 uM ADP	525	616	1031
2 uM ADP	524	579	896
Non	171	0	98

Figure 8

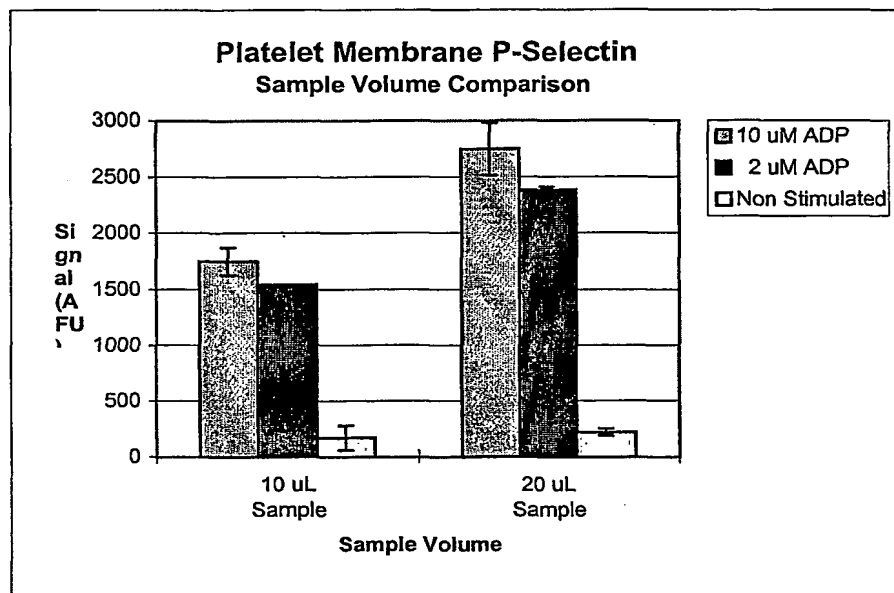

	<u>10 uL</u>	<u>20 uL</u>
10 uM ADP	1749	2747
2 uM ADP	1544	2390
Non	170	219

INTERNATIONAL SEARCH REPORT

Int'l Appl. No.
PCT/US 01/2132A. CLASSIFICATION OF SUBJECT MATTER
IPC 7 G01N33/86 G01N33/543

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 G01N

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

BIOSIS, EPO-Internal, WPI Data

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	DE 42 28 395 A (LOELIGER CHRISTOPH CORNELIUS D) 3 March 1994 (1994-03-03)	1-11, 13-16, 18-20
Y	the whole document	12, 17
X	WO 00 25140 A (ACCUMETRICS INC) 4 May 2000 (2000-05-04)	1-5, 10, 11, 13-16, 19 6-9, 17, 18
Y	abstract page 14, line 24 - line 26 page 16, line 23 - page 17, line 18 page 18, line 20 - line 26 claims 1-14; examples 4-7	
	----- -/--	



Further documents are listed in the continuation of box C.



Patent family members are listed in annex.

* Special categories of cited documents:

- *A* document defining the general state of the art which is not considered to be of particular relevance
- *E* earlier document but published on or after the international filing date
- *L* document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
- *O* document referring to an oral disclosure, use, exhibition or other means
- *P* document published prior to the international filing date but later than the priority date claimed

- *T* later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
- *X* document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
- *Y* document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.
- *&* document member of the same patent family

Date of the actual completion of the international search

20 November 2001

Date of mailing of the international search report

17/12/2001

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax (+31-70) 340-3016

Authorized officer

Stricker, J-E

INTERNATIONAL SEARCH REPORT

Int'l Application No
PCT/US 01/24132

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category *	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	WO 90 04019 A (BAXTER INT) 19 April 1990 (1990-04-19) claims 22-32, 36-39 ---	6-9, 18
Y	JENSEN MORTEN KROGH ET AL: "Increased platelet activation and abnormal membrane glycoprotein content and redistribution in myeloproliferative disorders." BRITISH JOURNAL OF HAEMATOLOGY, vol. 110, no. 1, July 2000 (2000-07), pages 116-124, XP002183412 ISSN: 0007-1048 abstract page 118, column 1, paragraph 3 -column 2, paragraph 1 ---	12
Y	DATABASE BIOSIS 'Online! BIOSCIENCES INFORMATION SERVICE, PHILADELPHIA, PA, US; 1998 CHEN RIYAN ET AL: "Cytoskeletal changes in platelets induced by thrombin and phorbol myristate acetate (PMA)." Database accession no. PREV199900308212 XP002183414 abstract & CELL BIOLOGY INTERNATIONAL, vol. 22, no. 6, 1998, pages 429-435, ISSN: 1065-6995 ---	17
Y	DATABASE BIOSIS 'Online! BIOSCIENCES INFORMATION SERVICE, PHILADELPHIA, PA, US; 1996 MUSTONEN PIRJO ET AL: "Epinephrine augments platelet recruitment to immobilized collagen in flowing blood: Evidence for a von Willebrand factor-mediated mechanism." Database accession no. PREV199698722467 XP002183415 abstract & THROMBOSIS AND HAEMOSTASIS, vol. 75, no. 1, 1996, pages 175-181, ISSN: 0340-6245 ---	17
A	WHITE JAMES G ET AL: "Glycoprotein Ib is homogeneously distributed on external and internal membranes of resting platelets." AMERICAN JOURNAL OF PATHOLOGY, vol. 155, no. 6, December 1999 (1999-12), pages 2127-2134, XP002183413 ISSN: 0002-9440 cited in the application abstract page 2128, column 2 -----	

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No.

PCT/US 01/21132

Patent document cited in search report		Publication date	Patent family member(s)	Publication date
DE 4228395	A	03-03-1994	DE 4228395 A1	03-03-1994
WO 0025140	A	04-05-2000	AU 1218900 A	15-05-2000
			EP 1123510 A1	16-08-2001
			WO 0025140 A1	04-05-2000
WO 9004019	A	19-04-1990	AU 4492789 A	01-05-1990
			CA 1335181 A1	11-04-1995
			DE 68917736 D1	29-09-1994
			DE 68917736 T2	04-05-1995
			EP 0438520 A1	31-07-1991
			JP 7057177 B	21-06-1995
			JP 4500008 T	09-01-1992
			WO 9004019 A1	19-04-1990
			US 5536475 A	16-07-1996

THIS PAGE BLANK (USPTO)

**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☐ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.

THIS PAGE BLANK (USPTO)